


Physiological Signal Record (NCCU)

Step 1: After setting up with pulse oximetry and respiration belt, you should check waveform of respiration and pulse.



Step 2: Open Windows Command shell window on the MRI console



→ Type: *ideacmdtool*

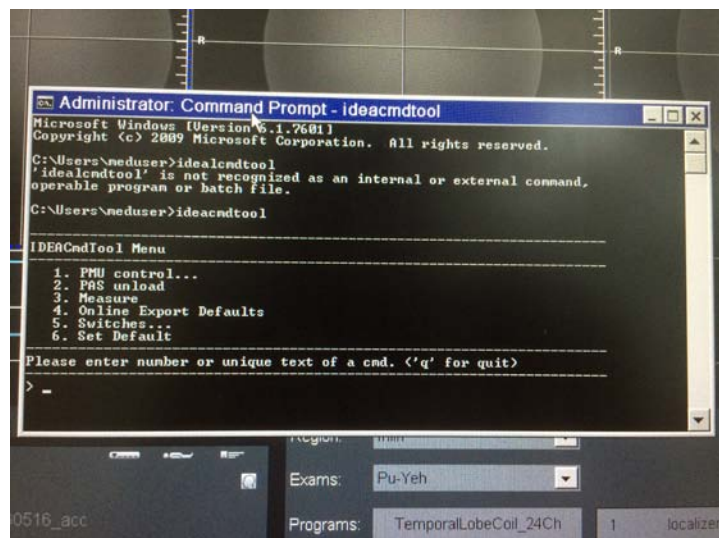

Step 3: Open PMU logging tool

→ Choose: 1 (PMU control)

→ Choose: 2 (PMU Signal logging)

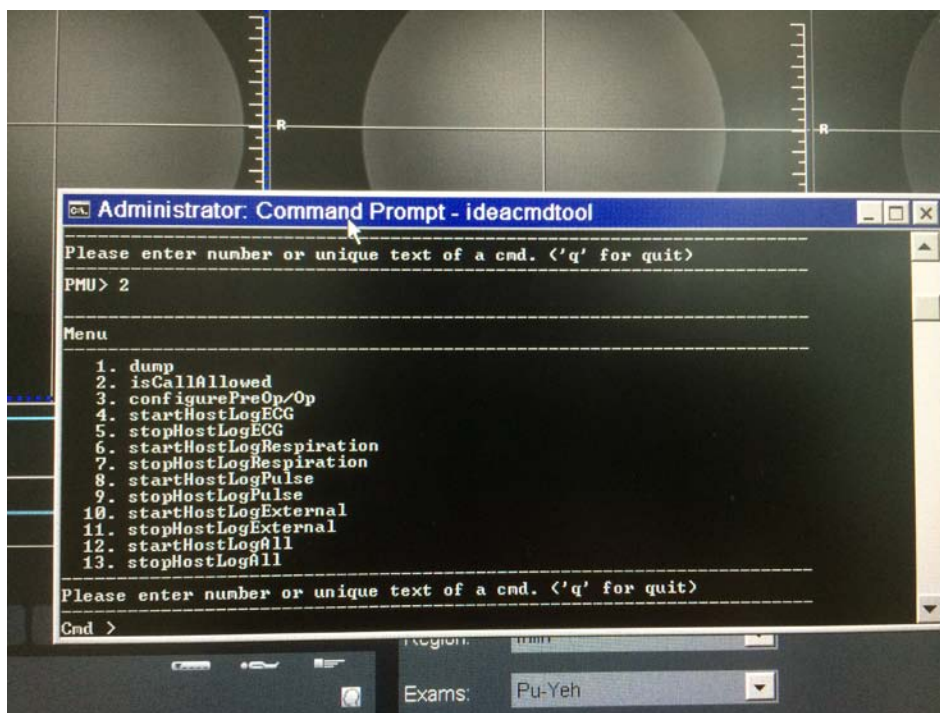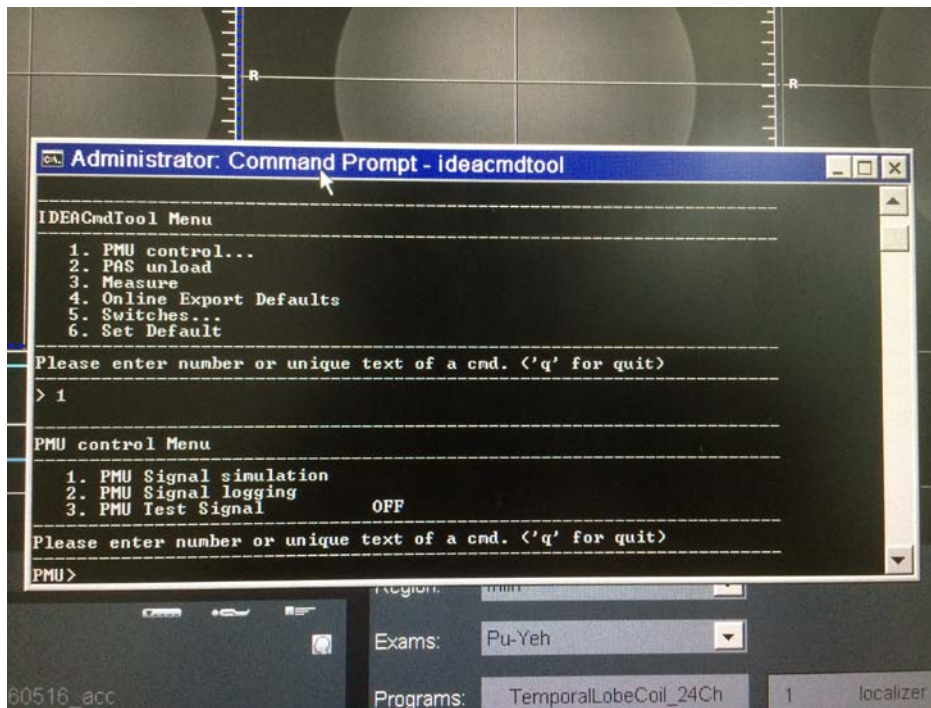



Step 4: Start the signal logging

→ Choose: **12** (startHostLogAll)

→ Enter log file name (without path and extension)

Ex: 161024_XXX_subjX_runX, and press “enter” to start recording

→ When the scan is completed → **enter: 13**

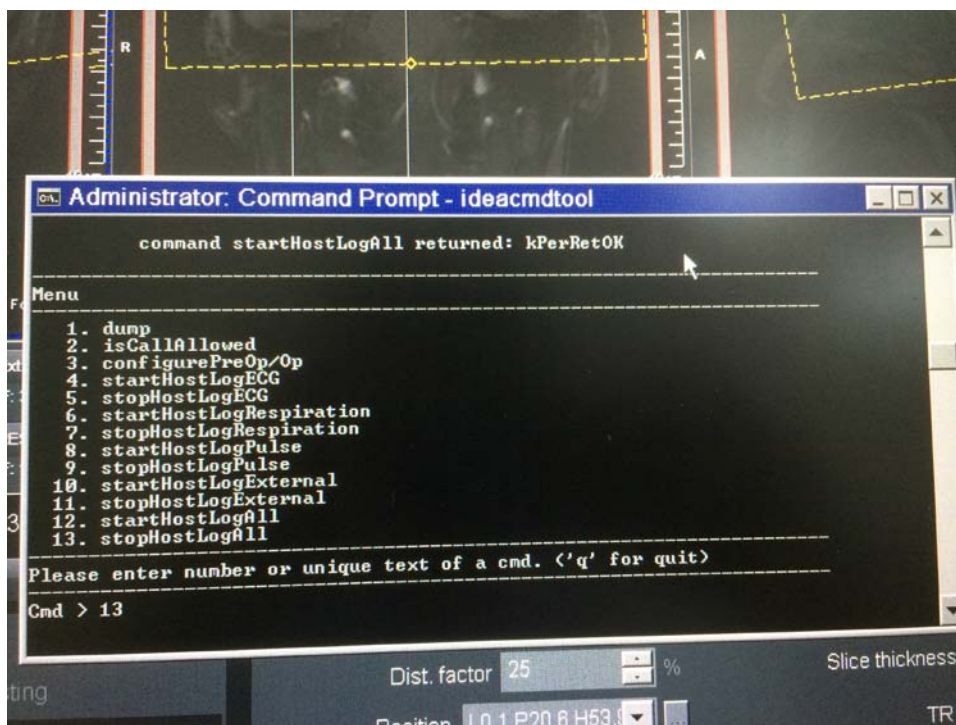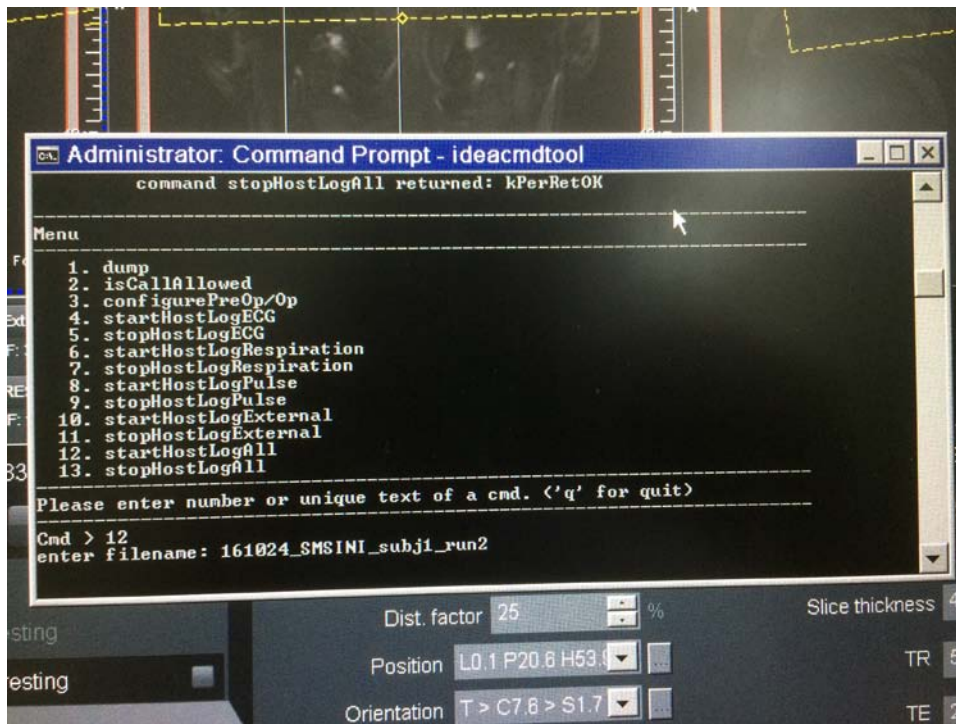

Step 5: The resulting log files can be found in the directory **c:\MedCom\log**

161024_XXX_subjX_runX.ext + 161024_XXX_subjX_runX.resp

+ 161024_XXX_subjX_runX.puls + 161024_XXX_subjX_runX.ecg

