


## Physiological Signal Record (Trio VB 17)

Step 1: After setting up with pulse oximetry and respiration belt, you should check waveform of respiration and pulse.



Step2: Open Windows Command shell window on the MRI console → Choose Menu item “Run...” → Type: *telnet MPCU*


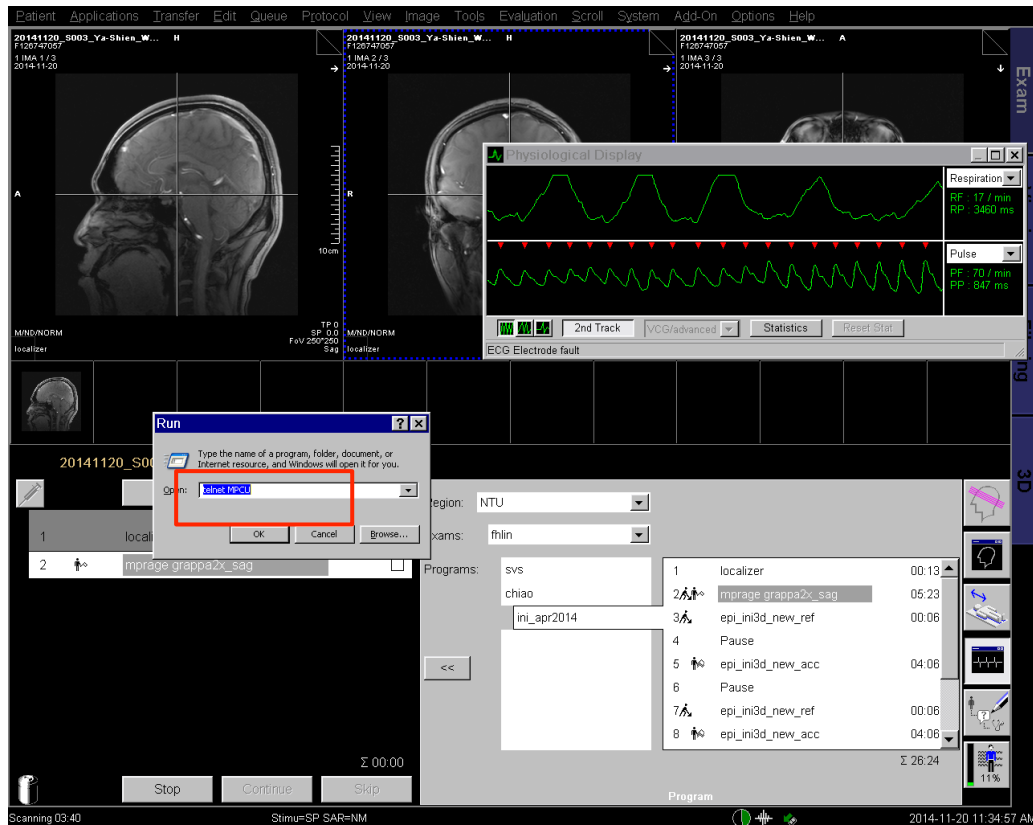


Step3. Start the signal logging:

Type: **fMenu** → enter 3 (choose + 3 - PMU ...)

→ enter 3 (choose +3 – Signal Logging ...)

→ enter 9 (choose +9 – Start all signals ...)

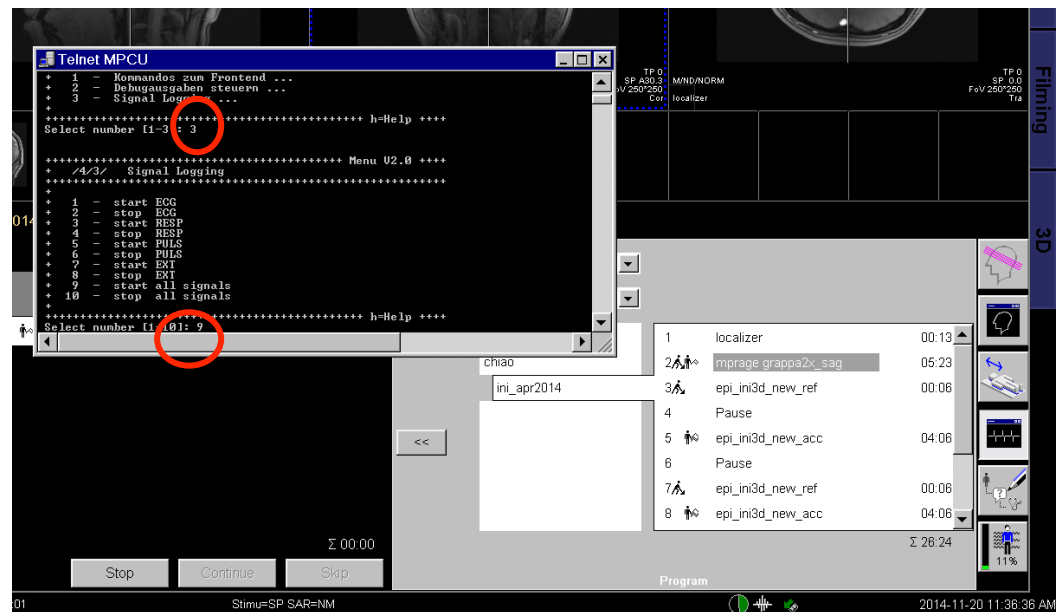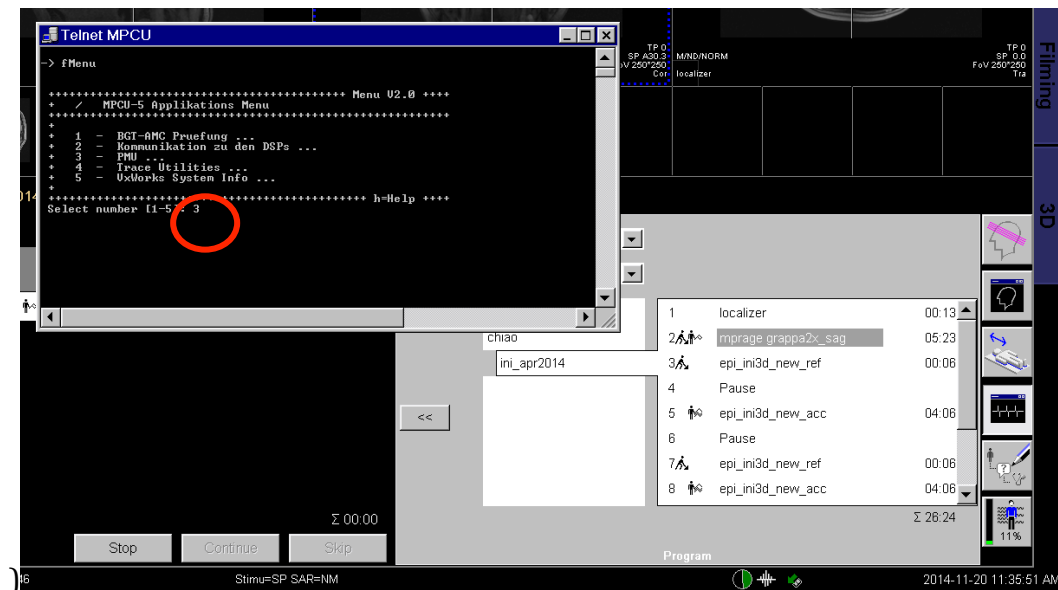

Step 4. Enter log file name (without path and extension):

Ex: 20141120\_XXX\_1(1 means the first run ), and press “enter” to record

❖Record PMU first →Run your MRI pulse sequence

When the scan is completed → **enter: 10**(choose +10– Stop all signals ...)



Step 5. The resulting logfiles can be found in the directory **c:\MedCom\log**

20141120\_XXX\_1.ext + 20141120\_XXX\_1.resp + 20141120\_XXX\_1.puls + 20141120\_XXX\_1.ecg

